Cloud-Clone Corp.

APA014Hu02 100µg Active Bone Morphogenetic Protein 4 (BMP4) Organism Species: Homo sapiens (Human) Instruction manual

#### FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

1th Edition (Apr, 2016)

### [PROPERTIES]

**Source:** Prokaryotic expression.

Host: E. coli

Residues: Ser22~Arg408

Tags: N-terminal His-tag

**Purity: >98%** 

**Buffer Formulation:** 20mM Tris, 150mM NaCl, pH8.0, containing 0.05% sarcosyl and 5% trehalose.

Applications: Cell culture; Activity Assays; In vivo assays.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.0

Predicted Molecular Mass: 48.0kDa

Accurate Molecular Mass: 48kDa as determined by SDS-PAGE reducing conditions.

# [ <u>USAGE</u> ]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

# Coud-Clone Corp.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### [SEQUENCE]

SHASLIPET GKKKVAEIQG HAGGRRSGQS HELLRDFEAT LLQMFGLRRR PQPSKSAVIP DYMRDLYRLQ SGEEEEQIH STGLEYPERP ASRANTVRSF HHEEHLENIP GTSENSAFRF LFNLSSIPEN EVISSAELRL FREQVDQGPD WERGFHRINI YEVMKPPAEV VPGHLITRLL DTRLVHHNVT RWETFDVSPA VLRWTREKQP NYGLAIEVTH LHQTRTHQGQ HVRISRSLPQ GSGNWAQLRP LLVTFGHDGR GHALTRRRA KRSPKHHSQR ARKKNKNCRR HSLYVDFSDV GWNDWIVAPP GYQAFYCHGD CPFPLADHLN STNHAIVQTL VNSVNSSIPK ACCVPTELSA ISMLYLDEYD KVVLKNYQEM VVEGCGCR

# [ACTIVITY]

BMP4 (Bone morphogenetic protein 4) is a member of the bone morphogenetic protein family, which is involved in bone and cartilage development, specifically tooth and limb development and fracture repair. It has been proven that HJV (Hemojuvelin) acts as a coreceptor of BMPs, including BMP4; therefore, a binding ELISA assay was constructed to detect the association of recombinant human BMP4 with recombinant human HJV. Briefly, BMP4 were diluted serially in PBS, with 0.01%BSA (pH 7.4). Duplicate samples of 100uL were then transferred to HJV-coated microtiter wells and incubated for 2h at 37°C. Wells were washed with PBST and incubated for 1h with anti-BMP4 pAb, then aspirated and washed 3 times. After incubation with HRP labelled secondary antibody, wells were incubated 15-25 minutes at 37°C. Finally, add 50µL stop solution to the wells and read at 450nm immediately. The binding activity of BMP4 with HJV was shown in Figure 1 and this effect was in a dose dependent manner.

Coud-Clone Corp.

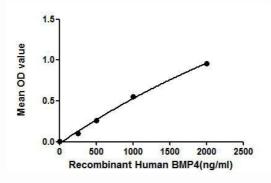


Figure 1. The binding activity of BMP4 with HJV.

### [IDENTIFICATION]


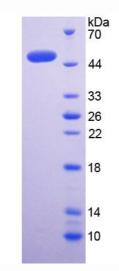


Figure 2. SDS-PAGE

Sample: Active recombinant BMP4, Human



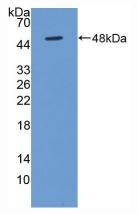

Sample: Recombinant BMP4, Human;

Antibody: Rabbit Anti-Human BMP4 Ab (PAA014Hu02)