## Cloud-Clone Corp.

PAA936Mu01 Polyclonal Antibody To Phosphoenolpyruvate Carboxykinase 1, Soluble (PCK1) Organism Species: Mus musculus (Mouse) Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

#### [PRODUCT INFORMATION]

Immunogen: PCK1, Mouse Clonality: Polyclonal Host: Rabbit Immunoglobulin Type: IgG Purification: Affinity Chromatography. Applications: WB, ICC, IHC-P, IHC-F, ELISA Concentration: 200µg/mL UOM: 100µg

#### [IMMUNOGEN INFORMATION]

Immunogen: Recombinant PCK1 (Ile312~Ala563) with two N-terminal Tags,

His-tag and T7-tag expressed in E.coli.

Accession No.: RPA936Mu01

#### [ ANTIBODY SPECIFITY ]

The antibody is a rabbit polyclonal antibody raised against PCK1. It has been selected for its ability to recognize PCK1 in immunohistochemical staining and western blotting.

## [APPLICATIONS]

Western blotting: 1:50-400 Immunocytochemistry in formalin fixed cells: 1:50-500 Immunohistochemistry in formalin fixed frozen section: 1:50-500 Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

# Cloud-Clone Corp.

### [ CONTENTS ]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

#### [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant PCK1 (Ile312~Ala563) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.


5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.


**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.

## [ <u>STORAGE</u> ]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.



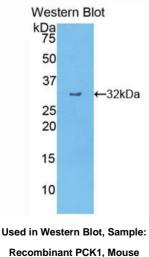

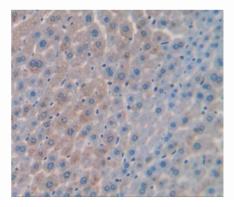

Used in DAB staining on fromalin fixed paraffin- embedded liver tissue