

PAD665Hu01 Polyclonal Antibody to Acetyl Coenzyme A Acetyltransferase 2 (ACAT2) Organism Species: Homo sapiens (Human) Instruction manual

> FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

> > 12th Edition (Revised in Aug, 2016)

1304 Langham Creek Dr. Suite 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-cloue.us | mail@cloud-clone.as Export Processing Zone, Wuhan, Hubei 436056, PRC | 0086-000-880-0687 | www.cloud-clone.com | mail@cloud-clone.com

# Cloud-Clone Corp.

## [PROPERTIES]

Source: Polyclonal antibody preparation Host: Rabbit Purification: Antigen-specific Affinity Chromatography. Traits: Liquid Concentration: 200µg/mL UOM: 100µg Applications: WB; ICC; IHC-P; IHC-F; ELISA; IP; IF; FCM.

### [IMMUNOGEN]

Immunogen: Recombinant ACAT2 (Met1~Glu397) expressed in *E.coli*. Accession No.: RPD665Hu01

### [APPLICATIONS]

Western blotting: 0.5-2ug/ml Immunocytochemistry in formalin fixed cells: 5-20ug/ml Immunohistochemistry in formalin fixed frozen section: 5-20ug/ml Immunohistochemistry in paraffin section: 5-20ug/ml Enzyme-linked Immunosorbent Assay: 0.05-2ug/ml Optimal working dilutions must be determined by end user.

## [FORMULATION]

**Form & Buffer:** Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

## [QUALITY CONTROL]

**Content:** The quality control contains recombinant ACAT2 disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

# Coud-Clone Corp.

**Loading Buffer:** 100mM Tris(pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB, 50mM DTT and 0.02% NaN<sub>3</sub>.

#### [STORAGE AND STABILITY]


Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.


Aliquot and store at -20°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### [IDENTIFICATION]

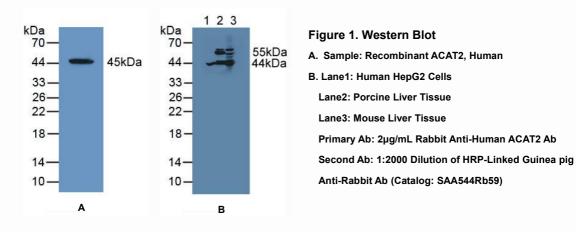

# Coud-Clone Corp.

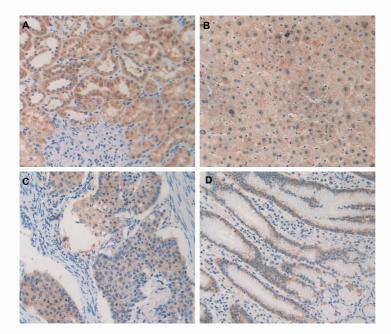

# Figure 2. DAB staining on IHC-P

Samples:

- A. Human Kidney Tissue
- B. Human Liver Tissue
- C. Human Lung Cancer Tissue
- D. Human Stomach Tissue

1304 Langham Creek Dr. Soire 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-cloue.ox | mail@cloud-clone.ax Export Processing Zone, Wuhan, Hubei 430056, PRC | 0086-000-880-0687 | www.cloud-clone.com | mail@cloud-clone.com