PAH339Hu01 Polyclonal Antibody to Periostin (POSTN) Organism Species: Homo sapiens (Human) *Instruction manual* 

# FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

1304 Langham Creek Dr. Soite 226, Houston, TX 77094, USA | 001-888-960-7402 | www.cloud-cloue.ox | mail@cloud-clone.ax Export Processing Zone, Wuhan, Hubei 436056, PRC | 0086-000-880-0687 | www.cloud-clone.com | mail@cloud-clone.com

# Cloud-Clone Corp.

#### [PROPERTIES]

Source: Polyclonal antibody preparation Host: Rabbit Purification: Antigen-specific Affinity Chromatography. Traits: Liquid Concentration: 200µg/mL UOM: 100µg Applications: WB; ICC; IHC-P; IHC-F; ELISA; IP; IF; FCM.

#### [IMMUNOGEN]

Immunogen: Recombinant POSTN (Pro97~Leu230) expressed in *E.coli*. Accession No.: RPH339Hu01

#### [APPLICATIONS]

Western blotting: 0.5-2ug/ml Immunocytochemistry in formalin fixed cells: 5-20ug/ml Immunohistochemistry in formalin fixed frozen section: 5-20ug/ml Immunohistochemistry in paraffin section: 5-20ug/ml Enzyme-linked Immunosorbent Assay: 0.05-2ug/ml Optimal working dilutions must be determined by end user.

#### [FORMULATION]

**Form & Buffer:** Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

#### [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant POSTN disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.

5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

### Coud-Clone Corp.

**Loading Buffer:** 100mM Tris(pH6.8), 1% SDS, 150mM NaCl, 50% glycerol, 0.02% BPB, 50mM DTT and 0.02% NaN<sub>3</sub>.

#### [STORAGE AND STABILITY]


Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

#### [ IDENTIFICATION ]

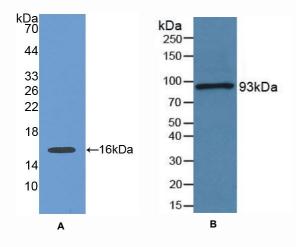


#### Figure 1. Western Blot A. Sample: Recombinant POSTN, Human

B. Human Serum

Primary Ab: 2µg/mL Rabbit Anti-Human POSTN Ab Second Ab: 1:2000 Dilution of HRP-Linked Guinea pig Anti-Rabbit Ab (Catalog: SAA544Rb59)

## Coud-Clone Corp.

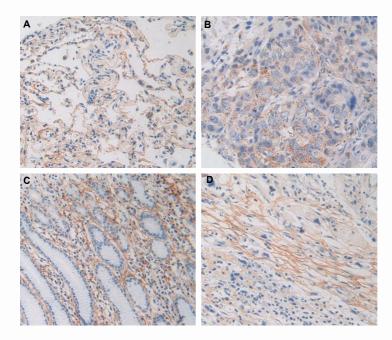

# Figure 2. DAB staining on IHC-P

Samples:

- A. Human Lung Cancer Tissue
- B. Human Breast Cancer Tissue
- C. Human Stomach Tissue
- D. Human Stomach Cancer Tissue

1304 Langham Creek Dr. Saite 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-cloue.us | mail@cloud-clone.as Export Processing Zone, Wuhan, Hubei 436056, PRC | 0086-000-880-0687 | www.cloud-clone.com | mail@cloud-clone.com