

PAH339Mu01

Polyclonal Antibody to Periostin (POSTN)

Organism Species: Mus musculus (Mouse)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

## [ PRODUCT INFORMATION ]

Immunogen: POSTN, Mouse Purification: Affinity Chromatography.

Clonality: Polyclonal Applications: WB, ICC, IHC-P, IHC-F, ELISA

Host: Rabbit Concentration: 200µg/mL

**Immunoglobulin Type**: lgG **UOM**: 100μg

## [ IMMUNOGEN INFORMATION ]

Immunogen: Recombinant POSTN (Ser500~Leu630) with N-terminal His-Tag

expressed in E.coli.

Accession No.: RPH339Mu01

## [ ANTIBODY SPECIFITY ]

The antibody is a rabbit polyclonal antibody raised against POSTN. It has been selected for its ability to recognize POSTN in immunohistochemical staining and western blotting.

## [ APPLICATIONS ]

Western blotting: 1:50-400

Immunocytochemistry in formalin fixed cells: 1:50-500

Immunohistochemistry in formalin fixed frozen section: 1:50-500

Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.



# [ CONTENTS ]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

# [ QUALITY CONTROL ]

**Content:** The quality control contains recombinant POSTN (Ser500~Leu630) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate. 5uL per well when used in enhanced chemilumescent (ECL).


**Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.


## [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.

## [IMAGES]

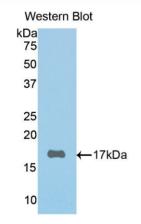

Used in Western Blot, Sample: Recombinant POSTN, Mouse

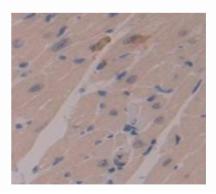

Used in DAB staining on fromalin fixed paraffin- embedded heart tissue